IT	Istruzioni d'uso.....	2
EN	Instructions for use.....	5
ES	Instrucciones de uso.....	8
DE	Gebrauchsanweisung.....	11
PT	Instruções de utilização.....	14
PL	Instrukcja użytkownika.....	17
FR	Instructions d'utilisation.....	20

USER MANUAL CWASH® EASY

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USER MANUAL CWASH® EASY

WARNING OF USE

The device must be used and maintained in an environment appropriate to the intended use and preferably stored at room temperature, not exposed to heat sources, to avoid possible general malfunctions, general structure deformation and/or damage to the electronic circuit.

The device is recommended and realized only and exclusively for the use according to which it was designed, therefore the manufacturer and the distributor are relieved of any type of liability for damage and/or problems caused by the incorrect use of the same.

Periodically check that the electrical cable is not damaged. Otherwise, please contact the assistance.

Do not use damaged or broken bases.

It is recommended to keep the device out of the reach of children. Not recommended for use by children under 3 years. For children under 6 years and for persons with reduced physical, sensory and mental capacity, it is recommended to use in the presence of an adult.

Do not ingest or chew.

In case of product fall, replace the bite before use even if there is no visible damage.

Do not place or store the charging base at a point where it could fall or be pushed into the bathtub or sink.

Do not immerse the device in water or other liquid.

Do not try to grab the device connected to the current in case of a fall in the water: immediately disconnect the plug from the socket.

Do not modify or repair the product, this may cause fire, electric shock or other injury.

Do not touch the power plug with wet hands: danger of electric shock.

Consult your dealer or contact service for repairs.

Do not disassemble the product.

Do not insert objects into the openings of the device.

When removing the power plug, pull the plug and not the cord.

Use the product only for the use described in this manual.

Do not use connections other than those recommended by the manufacturer.

If you are following a dental treatment, before using the device consult your dentist.

For the protection of the environment, do not throw the product away with normal household waste, but go to a disposal or recycling center in accordance with local regulations.

Do not dispose of in the environment after use.


INTRODUCTION

Cwash® Easy is the innovative dental cleaning device that doesn't need water and toothpaste. Small and easy to carry, Cwash® Easy can be used anywhere and anytime: will make you spend only 30 seconds of your time!

CONTENTS OF THE PACKAGE

The package contains:
1x Device Cwash® Easy
1x Gradual release bite
1x Closeable sack for bite
1x User manual

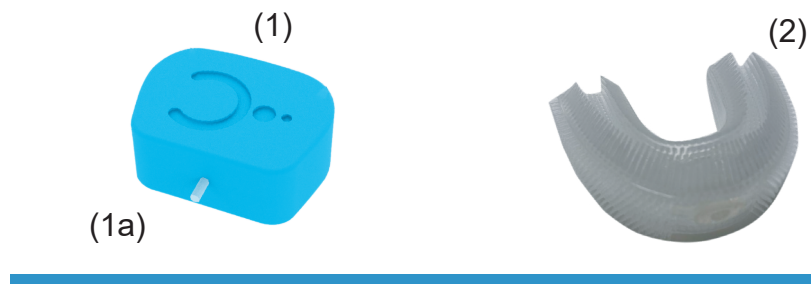ASSEMBLY AND INITIAL CONFIGURATION

Cwash® Easy is unassembled.

To proceed with its assembly it is necessary:

- Match the main device (1) with the gradual release bite (2) using the quick-release pin (1a);
- Apply slight pressure to facilitate connection between the two parts (1) and (2);
- At this point the device is ready to be used by clicking the "ON/OFF" button;
- The device will automatically shut off once the treatment is over or at the press of the "ON/OFF" button;
- After use, it is recommended to rinse the bite with water and store it in a dry and appropriate place.

CHANGING THE TREATMENT DURATION

Cwash® Easy is sold with the standard intensity.

To increase the treatment intensity, after having turned the device on:

- Press the "ON/OFF" button to increase the treatment intensity at medium level ;
- Press again the "ON/OFF" button to increase the treatment intensity at higher level.

Cwash® Easy at the end of the treatment, start again from the standard intensity.

REPLACEMENT OF BITE

After 100 cycles or 2 months of use, it is recommended to substitute the bite.

COMPLAINTS AND ASSISTANCE (FAQ)

For any complaints please contact your dealer or contact customercare@plusbiomedicals.com indicating the SERIAL NUMBER of the device which can be found inside the box cover.

To purchase replacement part, contact your dealer..

FAQs available on www.cwash-dental.com

PRIVACY AND PERSONAL DATA

Data processing is regulated with reference to the code on the protection of personal data and the GDPR (UE regulation 2016/679).

The full text is available on www.plusbiomedicals.com

Warranty of the product in compliance with the terms provided by Italian law.

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All rights reserved
Plus Biomedicals Srl